


Know and Maintain Your Folding Manual Wheelchair

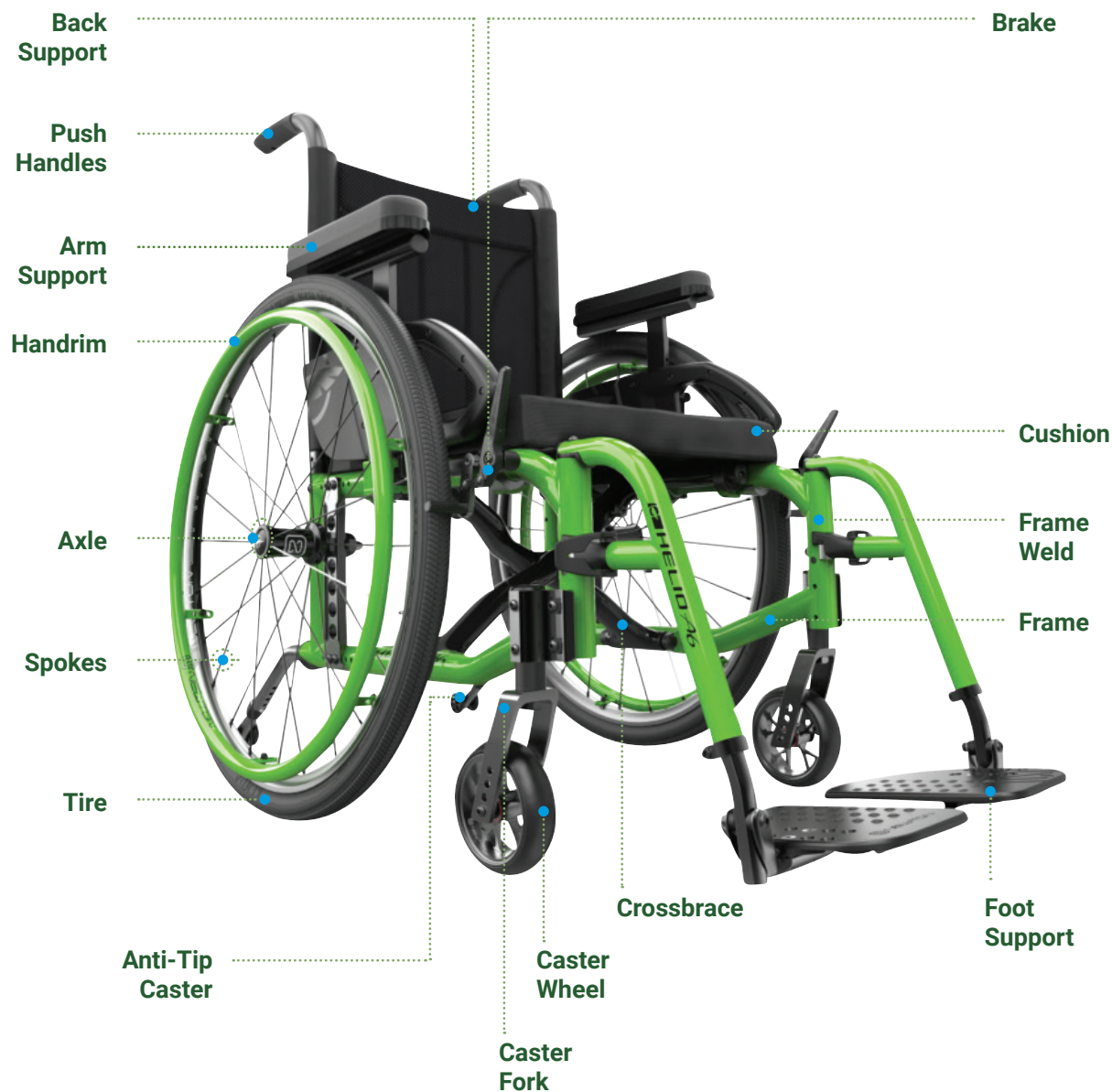

Maintenance Schedule

WEEKLY

Inspect cushion and cover

Inspect tire pressure (pneumatic ONLY)

MONTHLY

Casters

Clean and inspect caster bearings

Inspect anti-tip

Cushion

Wipe down cushion and cover

Frame

Inspect cross-brace, suspension, and welds

Wipe down frame

Rear Wheels

Inspect alignment, bearings, hand rim, brake, quick-release axle, spokes, and tires

Wipe down and lubricate rear wheel axle

Supports

Inspect arm, back, and foot supports

Inspect clothing guard

Inspect upholstery (including stitching & rivets)

Wipe down upholstery and supports